

DENTALAIRE™

2008 Catalog

Experience Excellence



800.866.6881 • www.dentalaireproducts.com

Meeting and Exceeding the Needs of
Veterinary Dentistry Since 1988

DENTALAIRE™

Veterinary Dental Division



Dentalaire has been providing our clients with quality products and great service since 1988. Our friendly and dedicated staff understands the needs of today's veterinary practice and are eager to help.

Our full line of products makes it easy to order and saves you time as well as money.

Good people, great service and quality products is what we're all about.

DENTALAIRE™

17150 Newhope Street, Suite 407 • Fountain Valley, CA 92708

1-800-866-6881

The Single Source for All Your Dental Needs.

Limited Warranty

DENTALAIRE warrants its products to be free from defects in materials and workmanship only. No other warranties are expressed or implied.

Air Compressors & Portable Dental Units:

All compressors and portable dental units including wall mount dental units are warranted for one (1) year. High Speed Turbines are warranted for six (6) months. High Speed & Low Speed Handpieces are warranted for one (1) year. Ultrasonic Scalers are warranted for two (2) years. Sonic Air Scaler Handpieces are warranted for one (1) year. Warranty is limited to repair, parts or replacement of item by DENTALAIRE. True Seal (Snap on and Screw Type) prophylaxis angles are warranted for one (1) year. **WARRANTY IS VOID IF ITEMS ARE CARELESSLY USED, IMPROPERLY MAINTAINED, IMPROPERLY INSTALLED, IMPROPERLY STERILIZED AND/OR MIS-USED. NO CLAIM FOR LABOR OR CONSEQUENTIAL DAMAGES WILL BE ALLOWED.**

Instruments:

All DENTALAIRE and Vetcare Dental Instruments are warranted for one (1) year. If any of our instruments break, chip, or bend during any dental procedure, Dentalaire will replace the instrument at no charge. Defective item must be returned to us for exchange. **WARRANTY IS VOID IF ITEMS ARE CARELESSLY USED, IMPROPERLY MAINTAINED, IMPROPERLY INSTALLED, IMPROPERLY STERILIZED AND/OR MIS-USED. NO CLAIM FOR LABOR OR CONSEQUENTIAL DAMAGES WILL BE ALLOWED.**

Terms:

Net 30 days on approved accounts. A service charge of 1-1/2% per month (order only) will be charged on past due accounts. All prices subject to change without notice. All prices are FOB Fountain Valley, California or Houston, Texas. 15% handling charge on all returned items, unless deemed defective. Returned items must be shipped prepaid. Items may not be returned without factory authorization. International orders must be accompanied by an irrevocable Letter of Credit. All shipping to be paid by customer.

When placing your order please include product number, description and quantity. Also please include your purchase order number (if applicable), shipping date required and preferred method of shipping.

Opened handpieces and equipment may not be returned for credit, but will be repaired or replaced in accordance with manufacturer warranties. Before opening handpieces or equipment, we suggest that you check the shipping container

Table of Contents

Curettes and Scalers.....	2-4	Burs and Diamond Discs.....	19-20
Instrument Trays/Instrument Sets	5	Endodontic Products	21
Elevators and Levitators.....	6-8	Q & A on Dental Stations.....	22
Woodward Extraction Kit.....	8	Equipment Comparison Chart.....	23
Extraction Forceps and Scissors.....	9	Dental Stations/Compressors.....	24-26
Rodent Dental Instruments.....	10	Ultrasonic Scalers	27-28
Ultrasonic Cleaner	10	MaxiSonic Combo and Micromotors	28
Prophy Products.....	11	Dental X-Ray, X-Ray Products.....	29-34
Mouth Gags and Retractors.....	12	Handpieces.....	35
Restorative Products	12-19	Sharpening System / Miscellaneous Items	36

HAND INSTRUMENTS FOR VETERINARY DENTISTRY

Understanding the purpose and use of hand instruments is basic to veterinary dental procedures. Each instrument is designed for a specific application from calculus removal to tooth extraction. The classification is determined by its use.

CURETTES are used to remove subgingival calculus, root planing and for curettage (removal of disease in the soft tissue of the periodontal pocket). They have one blade with two cutting edges. Both cutting edges may be used on the front and back of the tooth.



DTP11227 – Bamhart Curette 5/6



DTP11340 – McCall's 17/18 Curette



DTP11228 – Bamhart Curette 5/6 (Blue)



DTP11224 – Columbia 13/14

HAND INSTRUMENTS FOR VETERINARY DENTISTRY

PERIODONTAL PROBES AND EXPLORERS are detection instruments. The probe is used for measuring the depth of periodontal pockets. The explorer is used to detect calculus.



DTP11390 – Explorer/Probe 1/23



DTP00012 – Perlo-Probe, Single End 3-6-9-12



DTP00015 – Perlo-Probe, Single End 3-6-9-12-15

SCALERS are used to remove dental deposits from tooth surfaces above the gum line.



DTP11278 – Tartar Sickle/Chisel (Gold)



DTP11276 – Sickle Scaler H6/H7 (Green)



DTP11275 – Sickle Scaler H6/H7 (Economy)

HAND INSTRUMENTS FOR VETERINARY DENTISTRY

GRACEY CURETTES have only one cutting edge which is designed for deep subgingival calculus, root planing and for curettage of the periodontal pocket.



DTP11271 – Jacquette Feline Scaler 1/2 (Black)



DTP11274 – Morse Scaler 0/00 (Violet)



DTP11279 – Tartar Scraper Double End (Economy)
DTP11277 – Tartar Scraper Double End



DTP11226 – Gracey 7/8 Curette



DTP10545 – Dental Mouth Mirror Size # 3

SPECIAL "MINI" SERIES FOR SMALLER ANIMALS

SHORT SHANK HAND INSTRUMENTS



- DTP11282 M5/6 Gracey 5/6
- DTP11283 M2/3 Jacquette 2/3
- DTP11284 M0/00 Morse 0/00
- DTP11285 M6/7 Sickle Scaler H6/H7
- DTP11286 M5/6 Barnhart 5/6
- DTP11288 M15/16 Tartar Sickle/Chisel
- DTP11287 Complete Set of 6

These hand instruments are ideal for feline dentistry, as well as for scaling the teeth of rodents. The shorter 3/4" shank allows for easier access to the mouth and gives you much better stability. This in turn will greatly reduce hand fatigue.

**With Free
Instrument
Tray**

BASIC DENTAL INSTRUMENT SET-UP

DTP10612 – Includes all items listed (set of 6)

- Mouth Mirror
- Probe/Explorer
- Tartar Scraper
- Scaler H6/H7
- Gracey 7/8
- Barnhart 5/6



COLOR CODED INSTRUMENTS



- DTP11330 – Set of 6 w/Autoclavable Tray
- | | |
|----------------------|--------|
| Gracey Curette 5/6 | Silver |
| Barnhart Curette 5/6 | Blue |
| Morse Scaler 0/00 | Violet |
| Sickle Scaler H6/H7 | Green |
| Tartar Sickle/Chisel | Gold |
| Jacquette Scaler 1/2 | Black |

SHARPENING STONES



DTP10615 – Arkansas Sharpening Stone, Conical



DTP10375 – Arkansas Sharpening Stone, Flat



DTP10660 – Sharpening Stone Oil .5 oz.



DTP90605 – Acrylic Test Stick (3 pack)

INSTRUMENT TRAYS



DTP10925
Small Instrument
Tray Divided
Autoclavable



DTP10925L
Large Instrument
Tray (Flat)
Autoclavable



DTP10925D
Large Instrument
Tray (Divided)
Autoclavable

ELEVATORS

Straight Tip Elevators come with blades from 1mm-5mm wide which will cover all your needs except curved rooted teeth. SUPER SLIM tips are excellent for cats. MODIFIED TIPS are designed to guide it into the socket and cut the periodontal ligament. This prevents it from slipping out and causing damage to the gingival tissue. Makes elevation much easier.



Regular Tip



Modified Tip



Super Slim

DTP10331 – Feline Elevator



Super Slim

DTP10331SM – Feline w/V Cut



DTP20331 – Canine Elevator



DTP10330 – Std. 301 Elevator

DTP10226 –
Std. 304 Elevator

LEVITATORS

MODIFIED TIPS

The single ended feline and canine levitators have a "V" cut into the end of the instrument. This helps guide the tip into the socket, thus preventing it from slipping, causing damage to the gingiva.

CURVED ELEVATORS

Indispensable elevators for extracting curved-root teeth (i.e. canine teeth and incisors). The curved tips are designed to follow the concave and convex curvatures of the teeth helping to prevent breaking the roots.

DTP00403
Single Ended Canine
(Curved)DTP00027
Single Ended Canine
(Straight)DTP00204
Single Ended Feline
(Curved)DTP00028
Single Ended Feline
(Straight)

SERRATED ELEVATORS

These elevators come highly recommended by Dr. Peter Emilly. The serrations help to keep you from slipping or rolling off the root of the tooth. It does a great job in severing the periodontal ligament. It also helps grabbing the cementum of the tooth that is being sectioned.



Comes in three convenient sizes.

- DTP10338 – Feline Serrated Elevator – Small
- DTP10340 – Mid-Sized Serrated Elevator – Medium
- DTP20332 – Canine Serrated Elevator – Large

LUXATORS



These luxators feature a very thin edge and ergonomic bend, and are very useful when extracting medium to large size teeth.

1. DTP0730/1 – Straight 3mm
2. DTP0730/2 – Straight 5mm
3. DTP0730/3 – Curved 3mm
4. DTP0730/4 – Curved 5mm
5. DTP0730/5 – Luxators-Set of 4



- DTP10361 – Cryer Elevator #30
- DTP10362 – Cryer Elevator #31

These elevators are helpful in removing the caudal teeth in dogs.

WINGED ELEVATORS

"Winged Elevators" were designed and created by Dr. Robert Wiggs, DVM/AVDC. These elevators will help to make extractions easier in both dogs and cats. They are sharp and have a fine edge to them. All of the winged elevators have a unique spoon shape that hugs the tooth and gives the veterinarian some added leverage, by using it in a side to side twisting/rocking motion. They are also useful in between the teeth in the interdental space. These come in three styles: Original S/S (lightweight feel), Color Coded (aluminum, heavier feel) and Mini-Wings (for smaller hands).



- DTP10341 – #1 (1mm) Original
- DTP10342 – #2 (2mm) Original
- DTP10343 – #3 (3mm) Original
- DTP10344 – #4 (4mm) Original
- DTP10345 – Set of Four with case
- DTP10328 – #6 (6mm) Original
- DTP10329 – #8 (8mm) Original

- DTPRXLSET – Set of 4 USA Made Winged Elevators
- DTPRX9001L – #1 USA Made Winged Elevator
- DTPRX9002L – #2 USA Made Winged Elevator
- DTPRX9003L – #3 USA Made Winged Elevator
- DTPRX9004L – #4 USA Made Winged Elevator

- DTPRXSSET – Set of 4 USA Made Short Shank Winged Elevators
- DTPRX9001S – #1 USA Made Short Shank Winged Elevator
- DTPRX9002S – #2 USA Made Short Shank Winged Elevator
- DTPRX9003S – #3 USA Made Short Shank Winged Elevator
- DTPRX9004S – #4 USA Made Short Shank Winged Elevator



Winged Elevator #1
Used for feline teeth and root tips

Winged Elevator #2
Used for feline teeth and small dog roots

Winged Elevator #3
Used for cat canine and dog teeth with medium sized roots

Winged Elevator #4
Used for dog canine teeth and other larger sized roots



- DTP10356 – #1 (1mm) "Mini"
- DTP10357 – #2 (2mm) "Mini"
- DTP10358 – #3 (3mm) "Mini"
- DTP10359 – #4 (4mm) "Mini"
- DTP10360 – Set of Four with case



- DTP10347 – #1 (1mm) Color Coded Violet
- DTP10348 – #2 (2mm) Color Coded Gold
- DTP10349 – #3 (3mm) Color Coded Blue
- DTP10350 – #4 (4mm) Color Coded Black
- DTP10351 – Set of Four with case

ROOT TIP PICKS



Indispensable for removing small pieces of broken root tips.

- DTP10175 – Root Tip Pick #78 (Straight)
- DTP10176 – Root Tip Pick #79 (Left Angle)
- DTP10177 – Root Tip Pick #80 (Right Angle)
- DTP10178 – Root Tip Pick (set of 3)

BONE CURETTES



Single Ended

- DTP10352 – Molt #1
- DTP10353 – Molt #2
- DTP10354 – Molt #4

PERIOSTEAL ELEVATORS

Used to elevate the gingiva and mucosa away from the diseased pocket after an incision has been made with a scalpel on one or both sides of the pocket.



DTP10333 – Molt #2/4 Perio Elevator Double-end
This is the periosteal elevator in the very popular Woodward Extraction Kit. One end is 2mm, and the other end is 4mm. This can be used in both dogs and cats. It is an excellent all purpose perio elevator used for doing gum flaps. Each end has a slight bend or curve to it.



DTP10337 – Freer #4 Small Perio Elevator



DTP453110 – Periosteotome Double Ended
One of our most popular periosteal elevators. It has small, very delicate ends and is excellent for use in doing gum flaps in small dogs and cats. Both ends are completely straight.



DTP10339 – Molt #9 Perio Elevator Double-end



DTP10335 – West #2 Perio Elevator Double-end

WOODWARD EXTRACTION KIT



Woodward Extraction Kit
DTP10613

This extraction kit includes all the basic dental instruments required for extractions in dogs and cats. It was designed by Dr. Tony Woodward, DVM/AVDC from Colorado Springs, CO. "This kit has the instrumentation for almost all your extraction needs."

The kit contains the following items: DTP10613

One (set of 4) S/S Winged Elevators. The set includes 1-4 mm. These are sharp and have a fine edge to them. They have a unique spoon shape that hugs the tooth and gives you some added leverage, by using it in a side to side twisting or rocking motion.

One 3mm Curved Luxator and one 5mm Curved Luxator. These elevators are useful in medium to large-sized teeth. Mostly used in dogs. These are sharp and have a fine edge to them and sometimes will fit where the winged elevators won't.

One pair of 7" Curved Serrated Dean Scissors. These are helpful in lengthening and trimming the gingival flap as well as cutting sutures in the mouth.

One Molt 2-4 Perio Elevator. This is a double-ended instrument used for elevating gingival flaps in dogs and cats.

One Small Rongeur Forceps, 4.5". Used as an extraction forceps and bone contouring tool.

One Straight Root Tip Pick. Useful for removing fractured root tips below the gumline.

EXTRACTION FORCEPS

Small extraction forceps are excellent for grasping small teeth (especially in cats). All have spring action for easy manipulation.



DTP10645 – Small Extracting Forceps (4-1/2")
Beak tips are flat and close completely. Perfect for small cats teeth.



DTP11305 – Tartar Remover Forceps - Small (4-1/2")
DTP11307 – Tartar Remover Forceps - Large (5-1/2")
Designed to remove heavy deposits of calculus from the teeth.



DTP10644 – Small Rongeur Forceps (4-1/2")
DTP20644 – Medium Rongeur Forceps (5-1/2")
Useful for gently grasping flat surfaces of large teeth (i.e. canine) for extraction.



DTP20645 – Large Extracting Forceps (6-1/2")
Useful for large canine teeth, premolars and molars.



DTP20379 – Fragment Forceps - Angular (5")
Indispensable for broken root fragments and some small delicate teeth, such as incisors of cats.



DTP13914 – (5-1/2") Extraction Forceps with right angle
Excellent for small incisors and caudal teeth.

SCISSORS

These are useful in lengthening and trimming the gingival flap.



Gum & Tissue
"S" Shape Serrated
Scissors (4-1/2")
DTP20954



Serrated / Curved
Dean Scissors (7")
DTP30960

RODENT DENTAL INSTRUMENTS



DTP454500 – Rodent Mouth Gag
An essential tool when examining the mouth of a rabbit.



DTP454550 – Small Rodent Cheek Dilator
DTP454555 – Large Rodent Cheek Dilator
Allow visualization of molars and protects cheeks for instruments. Also useful as a gag for guinea pigs and chinchillas.



DTP454800 – Molar Rasp
Used to remove tooth spurs without damaging the delicate root structure.



DTP454600 – Molar Extraction Forceps
Designed to grip without crushing. Most abscesses around the rabbit head have their origin in molar tooth roots.



DTP454650 – Rodent Molar Cutters
Used to trim spurs and hook angles. Unique shape enables the surgeon to see what is being cut.

RABBIT INCISOR LUXATOR

Malocclusion in rabbits is a serious problem. Extraction is the only permanent solution. The roots of the incisors in rabbits are large relative to the jaw. This luxator is designed to match the curvature of the roots of the incisors and is double ended for different size incisors.



DTP454700 – Rabbit Incisor Luxator
DTP454850 – Rodent Dental Kit (Set of 7 of the above instruments)



DTP02443 – Rodent Bur Guard (right)
Is a soft tissue protector used on low speed handpieces.

PROPHY PRODUCTS



DTP00730 – VETCARE
Propphy Paste
250 gr. Jar Fruit Flavor

DTP00729 – VETCARE
Propphy Paste
Mint Flavor

DTP00735 – Propphy
Paste Finger Ring Cup
Stainless Steel



ORALDENT CLEANSING SOLUTION

DA100CS

0.1% Chlorhexidine solution. 8 oz.



INDIVIDUAL PROPHY PAST CUPS

Package of 100

DTP00731



TRUE-SEAL

Propphy Angle
(Sealed & Autoclavable)

DA10955 – Snap-on

ECONOMY

Propphy Angle
(Autoclavable)

DA551PSN – Snap-on



TRUE-SEAL

Propphy Angle
(Sealed & Autoclavable)

DA10956 – Screw-on

ECONOMY

Propphy Angle
(Autoclavable)

DA552PSC – Screw-on



PROPHY CUPS

DTP00720 – Snap-on Propphy Cup, box of 144

DTP00670 – Screw-on Propphy Cup, box of 144

DTP00665 – Latch Type Propphy Cup,



CLOSE UP



New 'TWISTER' DISPOSABLE PROPHY ANGLE

Propphy Cup reciprocates back and forth.
No paste splatter! No heat!
No more tangles from catching hairs!

DTP10967 – Twister Propphy Angle with Cup, Bag of 100

DTP00731 – Individual Propphy Paste Cups, Box of 100

DTP10965 – Standard Disposable Propphy Angle with Cup, box of 144



FLUORIDE QUICK-GEL

(Stannous Fluoride .4%)

Helps prevent cavities by strengthening the enamel.
Also has an anti-bacterial effect. 4.3 oz. bottle

DTP10522 – Cinnamon

DTP10524 – Berry

DTP10526 – Mint



FLUORIDE CONCENTRATE

1.23% Acidulated & phosphorilated
DTP10525 – 16 oz. bottle



SAFETY GLASSES (Clear) DTP360005



SPLASH SHIELD

(Protection from paste splatter/bacteria)
DTP11234 – Pkg. of 3



DISCLOSING SOLUTION

Used to locate calculus or plaque on teeth. Good to spot check thoroughness of prophylaxis. Paint on with Q-tip or brush. Plaque Disclosing Solution

DTP10695 - 2 oz.

MOUTH GAGS & RETRACTORS



DTP42251 – Dental Mouth Props (Assort. 3 sizes, bag of 50). Red, yellow & blue.



DTP42254
Large Mouth Gag 7"
DTP42250
Canine Mouth Gag 5"
DTP42256
Feline Mouth Gag 3"

SENN RETRACTOR

DTP42260
Indispensable for gingiva retraction as well as cheek retraction.



PLASTIC CHEEK RETRACTORS

Ideal for pulling back the lip corners of dogs for access to the caudal part of the mouth.
DTP42257 – Small 1-1/2" wide
DTP32590 – Large 2" wide



RESTORATIVE DENTAL PRODUCTS

ALGINATE SPATULA

DTP10213



ALGINATE POWDER

DTP10125



FLEXI-BOWL

DTP150110 – Green, 700ML



IMPRESSION TRAY

DTP00966 – Set of 8



TWINZ™ VPS VINYL POLYSILOXANE IMPRESSION MATERIAL SYSTEM

Unsurpassed level of detail and accuracy for perfect fitting restorations:

- Thixotropic to avoid dripping or slumping
- Exceptionally hydrophilic for accurate, void-free impressions and pours
- Excellent tear strength at mouth removal and for multiple pours
- High dimensional stability, no rush to pour models
- No taste or odor to trigger a salivary response

Manufactured and packaged by BISCO, a history of quality.

DTP40030 – TWINZ VPS (Vinyl Polysiloxane Impression Material) 50cc. Heavy Body, Fast Track, Color - Aqua. 40 second oral work time.

DTP40050 – TWINZ VPS (Vinyl Polysiloxane Impression Material) 50cc. Light Body, Fast Track, Color - Green. 40 second oral work time.

DTP41100 – TWINZ VPS (Vinyl Polysiloxane Impression Material) 2 x 350gm jars. Heavy Body, Putty, Color - Blue. 1.5 minute room temp. work time.



DTP41040 – TWINZ VPS (Bite Registration Material) 50cc. Heavy Body, Color - Brown. 30 second oral work time.

DTP21010 – TWINZ VPS, Auto-Mix Bisco Dispenser.

To be used with DTP40030, DTP40050, DTP41040

RESTORATIVE DENTAL PRODUCTS



BEFORE



AFTER

Repairing fractured teeth has never been easier. To properly use this kit, you'll need a Curing Light (page 16) and a Latch type Contra Angle (page 35)

BONDED SEALANT KIT**DTP99950**

This kit has all you will need to start applying bonded sealants for exposed vital dentin and bonding composite restorations. The kit comes complete including instructions. It contains the following items: (treats 30-44 teeth)

1. **Etch 37% Phosphoric Acid** (5gm syringe), is used to etch dentin and enamel prior to placing restorative materials.
2. **One Step**, which is a universal light cure dental adhesive.
3. **Fortify**, is a "unfilled resin" used for coating areas of sealed dentin and composite restorations.
4. **Super Snap Mini Kit**, is designed for easy, fast, and safe contouring, finishing, and polishing of fractured teeth and composites.
5. **Disposable Mixing Wells** (48), these are ideal for handling and applying dental materials dispensed in small quantities.
6. **Applicators, Disposable** (Light Blue) (100), ideal for precise application of dental materials.
7. **(2) Dura-White Stones** for use in high speed handpiece.

DTP2901P**BISBLOCK**

(Oxalate Dentin Desensitizer) 6ML.



Is an oxalate dentin desensitizer that is designed to prevent intra-tubular fluid movement which causes dentinal sensitivity. This is used after the acid etch and then before the dental adhesive, for e.g. One Step by Bisco.

Benefits are:

- BisBlock gives you deep penetration into the dental tubules
- BisBlock gives you long-lasting relief
- BisBlock does not interfere with adhesive bonding
- BisBlock is safe, tissue friendly, and resists dislodgement



This is a great product to use in dogs that have excessive wear on their teeth. It can be used in conjunction with the very popular Bonded Sealant Kit.

An insight on using bonded sealants...

Smoothing a fracture site and applying bonded resin sealants serves several purposes. Smoothing the fracture site decreases plaque adherence and calculus accumulation, and decreases soft tissue irritation from the fractured surface. Plugging the dentin tubules with the bonded sealants decreases sensitivity and helps prevent infection of the pulp chamber. After the application of bonded sealants, the tooth is able to provide a permanent "patch" on the inside through the deposition of reparative dentin. Typically, this procedure is only done once per fracture site.

RESTORATIVE DENTAL PRODUCTS



SUPER-SNAP® POLISHING SYSTEM

DTP24490

Designed for easy, fast and safe contouring, finishing and polishing of all micro filled and hybrid composites. The ultra-thin disks are very flexible, bending under pressure without breaking and allowing easier access to interproximal areas. You will find the complete range of disks in the Super-Snap Rainbow Technique Kit. The kit contains 180 disks (100 standard – 80 mini), 40 polystrips, 2 Dura-White Stones (shapes CN1 and FL2), 1 Composite Fine, Mid-Point and 4 CA mandrels – all presented in a semicircular rainbow organizer tray.



SUPER SNAP MINI KIT

DTP24491

This kit is designed for easy, fast and safe contouring, finishing and polishing of fractured teeth and composites.



COMPOSITE FINISHING KIT

DTP10510

Kit includes two Dura-Green and two Dura-White CN1 and FL2 stones, 6 Composite Fine Points.



ALL-ETCH

All-Etch 10% (Phosphoric Acid) 5gm. syringe. This is a semi-gel that can be used on dentin and enamel at the same time to produce micro-retentive surface characteristics for successful bonding.

DTP5501E – All Etch 10% 5gm. Syringe

DTP5363E – All Etch 10% 70gm.



ETCH 37% - DTP5503E

(Phosphoric Acid) 5gm. Syringes used to etch dentin and enamel prior to placing restorative materials. (most popular)

ETCH 37% - DTP5638E

(Phosphoric Acid) 70gm. Same as above.

BLUE TIPS FOR ACID ETCH

DTP5503F (bags of 30)

RESTORATIVE PRODUCTS

ALL BOND 2 (6ml. x 2)



This is a universal dental adhesive, which is still unsurpassed in its ability to bond to dentin, enamel, new or old composites and amalgams. It has a dual-cured primer, which means it can be used with or without a light cure unit. It also comes with an antimicrobial etchant that etches and disinfects at the same time.

DTP42500 – ALL BOND 2 w/Etch, DE Resin, Prebond Kit
DTP2510A – ALL BOND 2 (Refill) 6 ml. x 2



MIXING WELLS



Disposable, packs of 48
DTP80496 – These are ideal for adhesives, opaquers and all materials dispensed in small quantities.



MIXING PADS



2-3/8" x 2-3/8", 5 pads per pack
DTP80485 – Great to use with Vitrebond.



ONE STEP 6ml



This is a universal light cure dental adhesive, which bonds to a wide variety of dental substrates. Very useful in sealing exposed dentin tubules in fractured teeth.

DTP2001P



ONE STEP PLUS 6ml



This is a "filled version" of the clinically proven One Step. It is designed to bond to composite dentin, enamel, cast materials, treated porcelain, and set amalgam. A very strong bonding agent for composite restoratives.

DTP1102P



FORTIFY 5ml



Composite Surface Sealant Refill
Includes instructions and MSDS

DTP9031F– Fortify



This is a "unfilled resin" used for sealing of composite restorations. It fills in microscopic defects and improves the longevity of composite restorations. It also helps to prevent micro leakage. This product is also very useful for sealing exposed dentin in fractures, when used in combination with the One Step (light cure) or All Bond (self cure) dental adhesives.

RESTORATIVE DENTAL PRODUCTS

LITEX CURING LIGHT



DTP00680

ONE YEAR WARRANTY

- 75 Watts
- 12 Volt
- Tungsten halogen lamp
- Internally mounted supercool fan
- Automatic Timer allows preselection of cure times from 10 to 60 seconds with audible tone
- Autoclavable Fiber Optic Probe
- Wall mount or free standing
- 9mm general purpose probe and eye protection shield

MEGA LIGHT TURBO LED CURING LIGHT



DTP00683

This light is very portable, never heats up (LED), is quiet and very easy to read. No bulbs are needed.

ONE YEAR WARRANTY

Technical Specifications:

- Power supply: 100V to 240V
- Weight: 80gm
- Light Source: Super Blue LED
- Spectrum Range: 450-490nm
- Light Guide: 8mm angled optical fibers
- Standard Accessories:
 - 1 light guide (8mm), 1 handpiece and 1 handpiece holder
 - 1 oval protector, 2 mini-protectors

MEGA LIGHT CORDLESS



DTP00685

ONE YEAR WARRANTY

No Strings attached with the new Mega Light Cordless. The latest LED Technology. Two connection options, battery or direct. Digital display and easy to operate. The Mega Light CL comes equipped with a built-in light meter (radiometer), this feature allows operator to verify light's efficiency. Preset curing times of 5, 10, 20, 30 and 40 seconds.

BRUSH TIPS / BRUSH HANDLES



DTP80223 – Brush Tips, pack of 1000. Ideal for applying adhesives, activators, sealants and other liquids.
DTP80211 – Brush Handle is ideal for reaching any tooth surface and holds the above brush tips.

APPLICATORS (100 per pack)



DTP80241

DTP80242



DTP80247

DTP80242 – Applicator (disposable). Fine / Silver. For applying adhesives in large areas.
DTP80241 – Applicator (disposable). Regular / Light Blue. For applying adhesives in large areas.
DTP80247 – Brush Applicator. Dark Blue. These Endo Brushes are used when applying adhesives to a prepared canal, prior to cementation of a post.

RESTORATIVE DENTAL PRODUCTS

VITREBOND GLASS IONOMER LINER / BASE


DTP24492



Vitrebond is a dual-cure (both light and self-cure) glass ionomer useful as a base/liner, in endodontic restorations, and in the restoration of caries at the gingival margin. Fluoride release aids in preventing the extension of carious lesions.

VISIBLE LIGHT CURE CALCIUM HYDROXIDE LINER

DTP10575

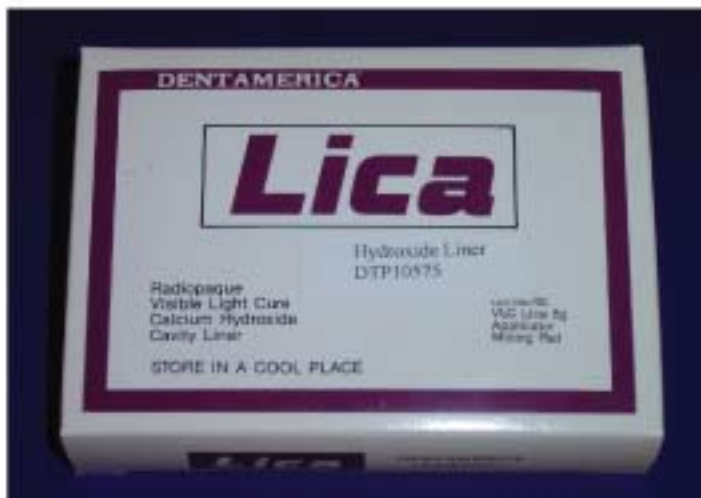

Used as a liner in cases of near pulp exposure, and in vital pulp therapy after placement of calcium hydroxide power on the exposed pulp tissues. Light cure polymerization is very convenient and shortens clinical procedure times.

GLASSIONOMER CEMENT

A self cure ionomer cement, useful as a base/liner and in the restoration of caries at the gingival margin. DTP10140

CEMENT SPATULA

DTP10212



AEHITEFLO COMPOSITE (LIGHT CURE / FLOWABLE)



This is easily injected into restoration sites. Benefits include radiopacity for easy identification, easily polished, very stain resistant, excellent wear resistance and strength. Used mainly for filling defects in non-occlusal areas.
 DTP671A2 – Shade A-2, refill 1.5gm. Syringe
 DTP671A3 – Shade A-3, refill 1.5gm. Syringe
 Light Cure / Flowable
 DTP673C2 – Shade C-2, refill 1.5gm. Syringe
 Light Cure / Flowable
 DTP674D3 – Shade D-3, refill 1.5gm. Syringe
 Light Cure / Flowable

RESTORATIVE DENTAL PRODUCTS

BIS-CORE Kit DTP1730B



2 syringe package, opaque. 1- Syringe dual cured catalyst 5gm.
1 syringe light cured base 5gm. Includes instructions and MSDS.

This core build up material is a "highly filled" composite available in two shades, (natural and opaque). It has a unique combination of viscosities allowing flexibility of use.



The base is light cured, the catalyst is a low viscosity self cured composite. It can be light cured to prevent slumping, while self cure finalizes the set.

DTP1731B – BIS-CORE REFILL
1 Syringe opaque catalyst 5gm.

DTP1732B – BIS-CORE REFILL
1 Syringe base 5gm.

DTP1734B – BIS-CORE REFILL
1 bottle liquid Catalyst 1.5gm.

WOODSON PFI #3 – Plastic Filling Instrument
DTP00210



BISFIL II KIT COMPOSITE DTP1210P



Self Cured Posterior Hybrid
Standard package - 1 Jar Base 7gm.
1 Jar Catalyst 7gm.
Accessories / Instructions / MSDS
Available in Cervical & Universal Shade

This is a self curing composite restorative 3 material that allows for restorative procedures without a light cure unit.

AELITE (LIGHT CURE / MICRO HYBRID) COMPOSITE



This is a very durable universal composite, which is ideal for all applications. Excellent for use in ALL CLASS of restorations. It offers you superior handling and resists slumping, which gives you greater control. It can be used as a core build up material or all alone by itself. This composite is polishable to an enamel like finish, while offering a high resistance to wear and stain. It can be used in high force areas of occlusion.

DTP814A2 - Aelite All Purpose Body - Composite (Shade A - 2)
DTP814A3 - Aelite All Purpose Body - Composite (Shade A - 3)
DTP814C2 - Aelite All Purpose Body - Composite (Shade C - 2)
DTP814D3 - Aelite All Purpose Body - Composite (Shade D - 3)

All 4 of these products come in a (.25 gm. syringe)



FLOUR PUMICE (FINE GRIT)
DTP10860 – 8 oz.

RESTORATIVE DENTAL PRODUCTS

**PRO-TEMP III
GARANT KIT**

Includes
dispensers and
cartridge with tips
DTP46201

**PRO-TEMP III
Replacement
Cartridge with tips
(21)**
DTP46246


This is a composite which is used to construct stabilizing occlusal splints for noninvasive mandibular fracture reduction. There are many uses such as: avulsion and luxation; symphyseal fracture reduction; turtle shell repair; bird beak fractures; bird beak alignment or occlusion correction. Also, to correct Base Narrow Lower Canines in dogs. Stronger, faster, easier to use than acrylics. Sets in 90 seconds.



CONSIL CERAMIC – Consil Bioglass, 2x2.5cc box
DTP10672

A unique bioactive ceramic for the regeneration of bone loss due to periodontal disease and tooth extraction. New bone production in 4 weeks with bending strength and stiffness in 3 months. Easy to use and mixes with patient's blood or sterile saline. Cost effective, you will get up to 10 procedures per box.

FRICTION GRIP DIAMOND BURS

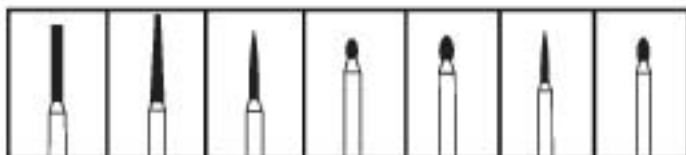

From Left to Right

- DTP7127C – Neodiamond Burs, Flat End Cylinder (Coarse)
- DTP1610C – Neodiamond Burs, Round End Taper (Coarse)
- DTP1610F – Neodiamond Burs, Round End Taper (Fine)
- DTP5128C – Neodiamond Burs, Flame (Coarse)
- DTP5128M – Neodiamond Burs, Flame (Medium)
- DTP5128F – Neodiamond Burs, Flame (Fine)
- DTP1908C – Neodiamond Burs, Egg (Coarse)
- DTP1908F – Neodiamond Burs, Egg (Fine)
- DTP1923M – Neodiamond Burs, Football (Medium)
- DTP128VF – Neodiamond Burs, Flame (Very Fine)
- DTP390VF – Neodiamond Burs, Egg (Very Fine)
- All Neodiamond Burs (5 per pack)

Diamond burs have a variety of uses, including crown preparation, shaping composite restorations, preparing sites for restorative materials, and smoothing bone in surgical sites prior to closure. One advantage is that they do not tend to catch and tear the soft tissues as do most other burs. The flat end cylinder burs are particularly useful for smoothing bone. A selection of grits allows for progressive refinement of the area you are working on. Diamond burs tend to have a much longer useful life than carbide burs. This selection represents the most popular shapes in veterinary dentistry. Total length never exceeds 22mm.

ADDITIONAL BURS

RA (Latch Type) BUR (5 per pack)

- DTP10036 – #700
- DTP10038 – #701
- DTP10042 – #702
- DTP10059 – #1/2
- DTP10060 – #2
- DTP10064 – #4
- DTP10066 – #6
- DTP10068 – #8
- DTP10040 – #33-1/2
- DTP10067 – #35
- Dentalaire Starter Bur Pack with
10 assorted RA (Latch type) Burs
- DTP20201 – Pk. of 10


HANDPIECE BURS (5 per pack)
(For Low Speed)

- | | |
|------------------|---------------|
| DTP10070 #700 | Handpiece Bur |
| DTP10071 #701 | Handpiece Bur |
| DTP10072 #702 | Handpiece Bur |
| DTP10073 #1/2 | Handpiece Bur |
| DTP10074 #2 | Handpiece Bur |
| DTP10075 #4 | Handpiece Bur |
| DTP10076 #6 | Handpiece Bur |
| DTP10077 #8 | Handpiece Bur |
| DTP10078 #33-1/2 | Handpiece Bur |
| DTP10079 #35 | Handpiece Bur |

FRICION GRIP BURS/ SURGICAL LENGTH BURS

Many types of burs are available. The cross-cut fissure bur (699, 700, 701, 701L, 702) is a good all-purpose bur that can be used to section teeth and gain endodontic access. It is also useful as a general purpose bur. The round bur (1/2, 1, 2, 3, 4, 6, 8) is used to gain endodontic access, to prepare cavities, and to mark the cutting depth in crown preparation. The inverted-cone bur (33-1/2, 34, 35, 37) is used for undercutting cavity preparation to create a mechanical interlock. The pear-shaped bur (330, 331, 332) is used for root canal access and for undercutting dentin.

ALL BURS COME IN PACKS OF 5

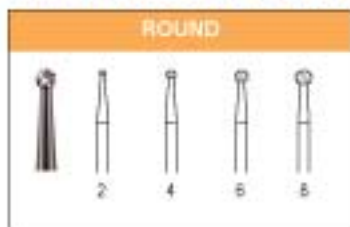

FRICION GRIP BURS (16mm.)

- DTP10046 – Round Bur #1/2
- DTP10047 – Round Bur #1
- DTP10048 – Round Bur #2
- DTP10050 – Round Bur #3
- DTP10052 – Round Bur #4
- DTP10054 – Round Bur #6
- DTP10056 – Round Bur #8
- DTP10201 – 10 Assorted Burs with Burr Block, includes #6, #8, #330, #332, #557, #558, #699, #700, #701, #702



- DTP10033 – Cone Shape #33-1/2
- DTP10034 – Cone Shape #34
- DTP10035 – Cone Shape #35
- DTP10037 – Cone Shape #37



- DTP10017 – Pear Shape #330
- DTP10018 – Pear Shape #331
- DTP10020 – Pear Shape #332
- DTP10021 – Cross Cut Fissure #557
- DTP10019 – Cross Cut Fissure #558
- DTP10025 – Cross Cut #699
- DTP10026 – Cross Cut #700
- DTP10028 – Cross Cut #701
- DTP10030 – Cross Cut #701L
- DTP10032 – Cross Cut #702

SURGICAL LENGTH BURS (25mm.)


- (Useful in big dogs with big teeth)*
- DTPFGSL2 – Surgical Length FG. #2
 - DTPFGSL4 – Surgical Length FG. #4
 - DTPFGSL557 – Surgical Length FG. #557
 - DTPFGSL558 – Surgical Length FG. #558
 - DTPFGSL6 – Surgical Length FG. #6
 - DTPFGSL701 – Surgical Length FG. #701
 - DTPFGSL702 – Surgical Length FG. #702
 - DTPFGSL8 – Surgical Length FG. #8



REPLACEMENT TIPS

FOR ULTRASONIC, MEGASONIC, AND ALL AIR SCALERS

TIP DIAGRAM

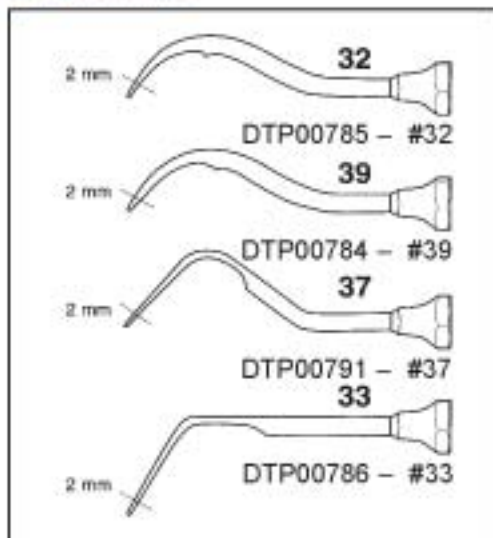

- DTP10880 – Cavitron Insert 25K (Fine)
- DTP10881 – Cavitron Insert 30K (Fine)
- DTP10879 – Cavitron Insert 25K (Wide)
- DTP100UN – Megasonic Tip (Universal)
- DTP100FL – Megasonic Tip (Flat)
- DTP100PR – Megasonic Tip (Perio)
- DTP100CT – Catheter Tip f/ Megasonic
- DTP100TW – Tip Wrench f/ Megasonic
- DTP10876 – Star Titan Universal Tip
- DTP10877 – Star Titan Perio Tip
- DTP10878 – Star Titan Rebuild Kit
- DTP10882 – Sonix, Lynx, Airstream (Universal)
- DTP10883 – Sonix, Lynx, Airstream (Perio)
- DTP10884 – Sonix, Lynx, Airstream (Sickle)
- DTP10870 – Tip for Air Scaler (Sickle)
- DTP10875 – Tip for Air Scaler (Universal)
- DTP00576 – Replacement Spindle for Air Scaler
- DTP10872 – Tip for Air Scaler (Perio)
- DTP10885 – Sonus V Stack with Nosecone
- DTP10886 – Universal Tip to Sonus V Stack
- DTP10887 – Flat Tip to Sonus V Stack
- DTP00784 – Universal Tip #39
- DTP00785 – Flat Tip #32
- DTP00786 – Special / Perio Tip #33
- DTP00791 – Interproximal #37

The #32 is for heavy duty scaling, the #39 is for fine or delicate scaling, the #33 and #37 are for subgingival scaling.

MAGNETIC BUR BLOCK ROUND



With cover - Holds 28 burs
DTP10204



DTP10205 Holds 7 burs

ENDODONTIC PRODUCTS

ENDODONTIC ORGANIZERS



BARBED BROACHES – (21mm)
Pack of 6 color coded. Used to remove root canal contents. DTP711330



Millimeter marked
ml. 029
Gutta Percha
DTP103301
240 pts. #15/80

Millimeter marked
as MMPP
Paper Points
DTP301203
240 pts. #15/80

HOLMSTROM PLUGGER/SPREADERS



- Special veterinary pluggers and spreaders designed by a veterinary dentist (Dr. Steve Holmstrom).
 - 35mm long, sizes #50, 65 & 90
 - 28mm long, sizes #20, 35
 - Especially useful for obstructions in larger canine teeth.
 - Available individually or in a set.
- DTP453120 – Holmstrom Plugger/Spreader #20
DTP453135 – Holmstrom Plugger/Spreader #35
DTP453150 – Holmstrom Plugger/Spreader #50
DTP453165 – Holmstrom Plugger/Spreader #65
DTP453190 – Holmstrom Plugger/Spreader #90
DTP453191 – Holmstrom Plugger/Spreader (Set of 5)

K-REAMERS



Set of 3 Assorted (60mm)
DTP01777 – #15-20-25
DTP01778 – #30-35-40
DTP01779 – #45-50-55
DTP01780 – #60-70-80

H-FILES



Set of 3 Assorted (60mm)
DTP01747 – #15-20-25
DTP01748 – #30-35-40
DTP01749 – #45-50-55
DTP01750 – #60-70-80

H-FILES



Set of 3 Assorted (120mm)
DTP01761 – #20-25-30
DTP01762 – #35-40-45
DTP01763 – #60-70-80

FINGER PLUGGERS



Set of 6 Assorted – (25mm)
DTP190662

- Used to vertically condense gutta percha in root canal procedures.
- 25 mm long, sizes #15-40.

FINGER SPREADERS



Set of 6 Assorted – (25mm)
DTP190676

- Used for laterally spreading gutta percha in root canal procedures.
- 25 mm long, sizes #15-40.
Color coded.

ENDO LOCKING PLIERS – DTP10570

Useful for handling paper point and gutta percha points.



60mm Paper Points (60)
DTP218291 – #15-20-25
DTP218292 – #30-35-40
DTP218293 – #45-50-55
DTP218294 – #60-70-80
DTP218295 – #90
DTP218296 – #100
DTP218297 – #110
DTP218298 – #120
DTP218299 – #130
DTP218300 – #140



60mm Gutta Percha (60)
DTP122203 – #15
DTP122204 – #20
DTP122205 – #25
DTP122206 – #30
DTP122207 – #35
DTP122208 – #40
DTP122209 – #45
DTP122210 – #50
DTP122211 – #55
DTP122212 – #60
DTP122213 – #70
DTP122214 – #80
DTP122215 – #90
DTP122216 – #100
DTP122217 – #110
DTP122218 – #120
DTP122219 – #130
DTP122220 – #140



DTP10230 – Zinc Oxide, 1.5 oz. (Powder)
DTP10153 – Eugenol (Liquid)

QUESTIONS & ANSWERS ON DENTALAIRE STATIONS

Q. How much horsepower do I need to run a Dentalaire unit?

A. 1/2 horsepower per dental station, which will run the high speed handpieces and low speed handpieces for hours of non-stop dental procedures.

Q. Do the dental stations come with high speed and low speed handpieces?

A. Yes, they do. We also include start-up supplies for polishing, a starter bur pack, and handpiece lubricant with all of our dental stations.

Q. Can Dentalaire customize a dental station to my certain needs or specifications?

A. Yes, we can. For example, we can add suction, a one horsepower motor, fiberoptics in the highspeed handpiece, ultrasonic scaling units, or almost anything else you desire. Since we are the manufacturer, we are happy to customize your unit to your desired individual specifications.

Q. What are the advantages of the Dentalaire compressor?

A. Dentalaire compressors are extremely quiet, odorless and are legendary for their workhorse construction.

Q. What is the warranty on the Dentalaire dental station?

A. They come with a one year warranty. The only exception is the piezo ultrasonic scaler, which comes with a two year warranty.

Q. Is Air Scaling the same as Piezo Ultrasonic Scaling?

A. No it is not. Air Scaling is "sonic" scaling rather than "ultrasonic", operating at approximately 6,000 cycles per second. Ultrasonic Scalers operate at up to 35,000 cycles per second. In general, the ultrasonic scalers are quieter and scratch the tooth surface less. Sonic scalers are better for subgingival use, as they generate less heat. Sonic scalers are air driven by the compressor, and place a large load on the compressor.

Q. What happens to me if my Dentalaire dental station breaks down 6 months after my warranty has expired?

A. Please call us at (800) 866-6881 to discuss your problem. Your satisfaction is our main goal. Many times we can trouble shoot over the phone to solve the problem.

Q. What needs to be done to keep the Dentalaire Dental Station running smoothly?

A. A. Always lubricate high speed handpieces and low speed handpieces before every use, except for the new DA700A Low Speed Handpiece, which is lube free. This takes only a few seconds.
 B. Use only distilled water in dental station to prevent build up in the lines and handpieces.
 C. Follow compressor maintenance guide chart. This takes only one minute per day.
 D. Make sure oil level is visible in the glass eye (at least half way up) on the back side of the compressor.
 Please note: Oil level is too high when you do not see any air bubble. If this happens, do not run compressor and call Dentalaire.

Q. What options do clinics usually add to their dental unit?

A. A. Push Button High Speed Handpiece.
 B. Fiberoptic Illumination in High-Speed Handpiece. This feature greatly improves your visibility and work quality.
 C. Adding a piezo scaler to any two handpiece unit
 D. Standard straight tubing or curled tubing that stretches.

Q. Does Dentalaire repair other companies equipment and handpieces?

A. Yes we do.

Q. Can I use other companies' handpieces in my new Dentalaire station?

A. Yes in most cases. Any standard "four-hole" handpiece will work with your Dentalaire station.

Q. My bur got stuck in my high speed handpiece. How can that happen?

A. By not putting the bur all the way in. This can cause small nicks in the shaft of the bur, causing it to become stuck. The tolerances for the shaft of the burs are very close, and even microscopic nicks in the shafts will not allow the burs to slide freely.

Q. How long does the light stay on in the high speed fiberoptic handpiece?

A. About 12 seconds after stopping the handpiece.

Q. Can I get suction in a Dentalaire system? Do I need to have suction?

A. Yes, you can get suction in a Dentalaire station. Most veterinarians work on a wet table, which eliminates the need for suction.

Q. When would I need a bigger compressor for my Dentalaire Dental Station?

A. If you desire suction or want to use a sonic scaling handpiece as your main method of scaling, you should consider a one horsepower silent air compressor. Suction and sonic scalers place a much larger load on a dental compressor.

EQUIPMENT COMPARISON CHART

Product #	Item	Hi Speed	Low Speed	Water & Air	Piezo Scaler	Air Scaler	Mobile Cart	Counter Top	Wall Mount	Table Mount	Fiber-optics	P.Button Bur	Compressor
DTP00494	Dentalaire II	X	X	X			X						
DTP00495	Dentalaire II	X	X	X			X						X
DTP00500	Dentalaire II	X	X	X				X					X
DTP00504	Dentalaire III	X	X	X	X			X					X
DTP00503	Dentalaire III with Cabinet	X	X	X	X			X					X
DTP00505	Dentalaire III	X	X	X		X		X					X
DTP00506	Dentalaire III	X	X	X		X		X			X		X
DTP00507	Dentalaire III	X	X	X		X	X						X
DTP00508	Dentalaire III	X	X	X		X	X						
DTP00509	Ultimate Vet Dental Station	X	X	X	X		X				X	X	X
DTP00510	Dentalaire III	X	X	X	X		X				X	X	
DTP00550	Cabinet to Countertop Unit												
DTP00600	Two Handpiece	X	X	X					X				
DTP00601	Two Handpiece	X	X	X					X				X
DTP00605	Three Handpiece	X	X	X		X			X				
DTP00606	Three Handpiece	X	X	X		X			X				X
DTP00607	Dentalaire Basic	X	X	X						X			
DTP00608	Dentalaire Basic	X	X	X						X			X

SCALERS / MICROMOTORS

DTP00772	Scalex "Classic" Ultrasonic Scaler: Has Cavitron Inserts, 25K
DTP00787	Scalex "Classic" Ultrasonic Scaler: Has Cavitron Inserts, 30K
DTP103RM	Ram 1000 Complete Prophy and Restorative Unit, 35,000 RPM Micromotor
DTP00780	Rotex 780 Micro Motor, 30,000 RPM
DTP00773	Pressurized Water Bottle System (Can be used on all Scalers)
DTP00776	Scalex 880 Scaler, Piezo Electric Scaler
DTP00770	Amdent US30 Scaler, Piezo Electric Scaler
DTP00789	MaxiSonic Scaler II Polisher Combo; Micromotor/Piezo Electric Scaler
DTP00766	Maxi-Sonic III w/ Reservoir
DTP00769	Swift Piezo Electric Scaler
DTP00765	Swift with Reservoir

All **Countertop Dental Systems** come with a **1000ml. water bottle**, which mounts on the inside of the compressor.

The **Mobile, Wall Mount, and Table Mount Systems** come with a quick connect **2000ml. water bottle**.

Distilled water is RECOMMENDED!

MOBILE CART DENTAL SYSTEMS



DTP00495 pictured

Dentalaire II Systems all have these basic features:

1. High Speed Handpiece
2. Low Speed Handpiece
3. Water and Air

There are six different models to choose from:

- DTP00494 – Dentalaire II without compressor
- DTP00495 – Dentalaire II with compressor
- DTP00507 – Dentalaire III with compressor and Air Scaler
- DTP00508 – Dentalaire III without compressor
- DTP00509 – The Ultimate Vet Dental Station (Dentalaire III) with compressor
See back cover
- DTP00510 – Dentalaire III without compressor
- DTP00777 – Upgrade to Piezo Scaler
- DTP0470 – Upgrade to Fiber Optic Push Button High Speed Handpiece.

Please refer to the equipment requirement chart on page 23 to see which features and options that certain dental systems have. Thank you!

DTP11310 – Compressor Oil 22 oz. (should be changed yearly.)

All Dental Stations include a FREE Dental Vet Package: 10 Friction Grip High Speed Burs, Stainless Steel Prophy Angle, Vetcare Prophy Paste (250gm jar), Prophy Cups (box of 144).

DENTALAIRE  1-800-866-6881




DTP00507 pictured

Dentalaire III Systems all have these basic features:

1. High Speed Handpiece
2. Low Speed Handpiece
3. Piezo OR Air Scaler
4. Water and Air

COUNTERTOP DENTAL SYSTEMS

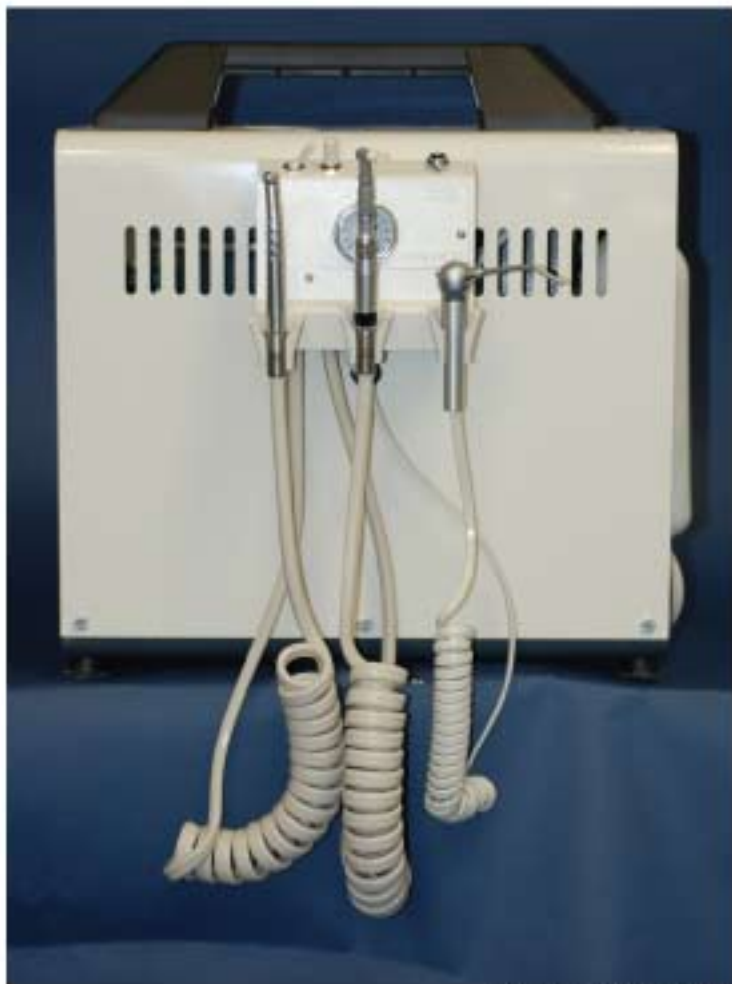

DTP00500 pictured

Dentalaire II Systems all have these basic features:

1. High speed Handpiece
2. Low speed Handpiece
3. Water and Air


There are five different models to choose from:

- DTP00500 – Dentalaire II with compressor
- DTP00504 – Dentalaire III with compressor and Piezo Scaler
- DTP00503 – Dentalaire III with compressor and mobile rolling cart
- DTP00505 – Dentalaire III with compressor and Air Scaler
- DTP00506 – Dentalaire III with compressor and fiber optics

Please refer to the equipment requirement chart on page 23 to see which features and options that certain dental systems have. Thank you!

DTP11310 – Compressor Oil 22 oz. (should be changed yearly.)

All Dental Stations include a FREE Dental Vet Package: 10 Friction Grip High Speed Burs, Stainless Steel Prophy Angle, Vetcare Prophy Paste (250gm jar), Prophy Cups (box of 144).

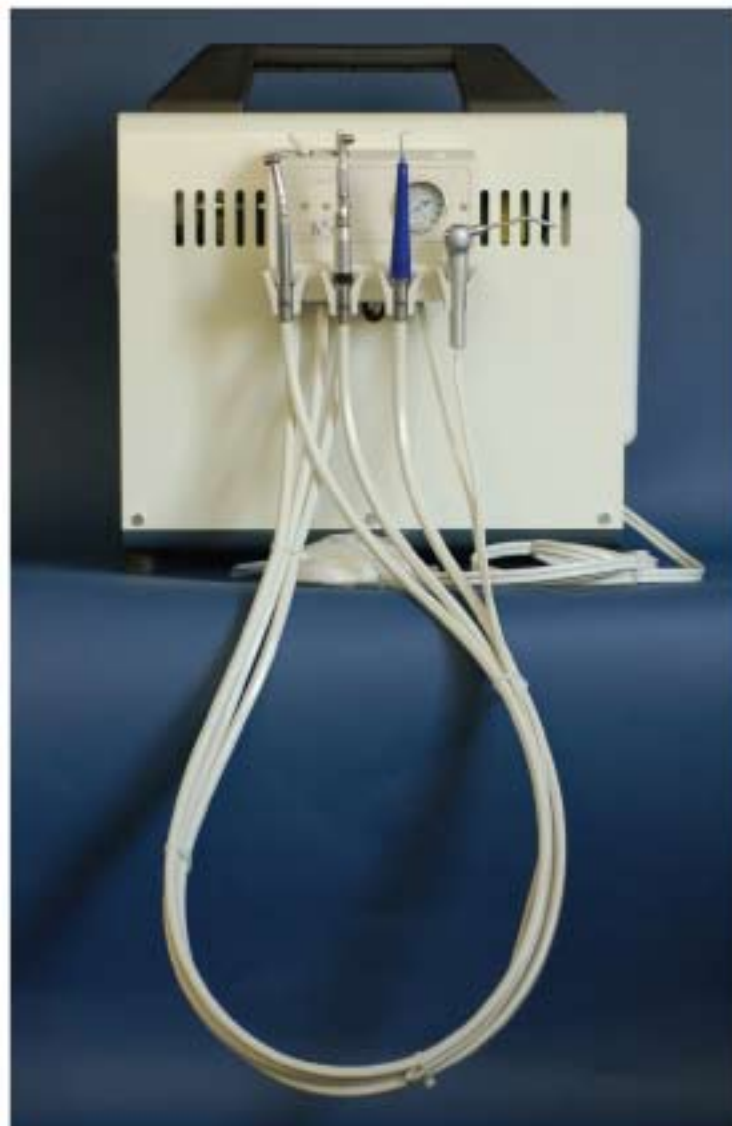

DTP00505 pictured

Dentalaire III Systems all have these basic features:

1. High speed Handpiece
2. Low speed Handpiece
3. Piezo OR Air Scaler
4. Water and Air

WALL MOUNT & TABLE MOUNT DENTAL SYSTEMS

WALL MOUNT DENTAL SYSTEMS



DTP00600

Two Handpiece System (pictured) has these basic features:

1. High speed Handpiece
2. Low speed Handpiece
3. Water and Air

Three Handpiece Systems all have these basic features:

1. High speed Handpiece
2. Low speed Handpiece
3. Air Scaler
4. Water and Air

DTP00777 – Upgrade Wall Mount System to Piezo Scaler

DTP00470 – Upgrade Wall Mount System to Fiber Optic Push Button High Speed Hand Piece

Wall mount systems available with 12, 14, 16 inch arms, and may be ordered as scissors (2) assembly, please specify when ordering.

DUAL USER SILENT AIR COMPRESSOR



Overall dimensions:
21.5" high, 27.5" wide, 13" deep
Tank: 6.5 Gallons
Capacity:
1 HP, 110V Motor delivers 4.3
CFM @ 80 PSI
110 V
Weight: 100 lbs

DA100-24

DENTALAIRE BASIC TABLE MOUNT DENTAL SYSTEMS



DTP00607

DTP00600 – Two Handpiece System without compressor

DTP00601 – Two Handpiece System with compressor

DTP00605 – Three Handpiece System without compressor, with Sonic Air Scaler

DTP00606 – Three Handpiece System with compressor and Sonic Air Scaler

DTP00607 – Dentalaire Basic without compressor, Two Handpiece System

DTP00608 – Dentalaire Basic with compressor, Two Handpiece System

SILENT AIR COMPRESSOR



Overall dimensions:
17" high, 7" deep, 16" wide
Tank: 2 Gallons
Capacity:
1/2 HP Motor delivers 2.1 CFM @ 80
PSI 110V
Weight: 72 lbs

DA50-9D



AVS

AVS Self Contained Air Powered Vacuum System. For optimal use, Dual User Silent air Compressor should be ordered, P/N DA100-24

All Dental Stations include a FREE Dental Vet Package: 10 Friction Grip High Speed Burs, Stainless Steel Prophy Angle, Vetcare Prophy Paste (250gm jar), Prophy Cups (box of 144).

DENTALAIRE SCALERS

SCALEX 880 PIEZO ULTRASONIC SCALER



DTP00776

Features:

- Self-tuning and load adjusting for maximum smooth operation
- Extremely gentle on patients
- Runs cool and quiet
- Compact size with slim lightweight handpiece.
- Includes 3 autoclavable tips, 2 universal and 1 flat.

Replacement Tips

- DTP00781 – Universal
- DTP00782 – Flat
- DTP00783 – Special/Perio

ONE YEAR WARRANTY

SCALEX-800 CAVITRON ULTRASONIC SCALER



DTP00772

Available in 25K or 30K
Complete with Cavitron Insert

DTP00772 25K

DTP00787 30K

Replacement Cavitron Inserts

- DTP10880 – 25K (Fine)
- DTP10881 – 30K (Fine)
- DTP10879 – 25K (Wide)

ONE YEAR WARRANTY

AMDENT US30 PIEZO ULTRASONIC SCALER



DTP00770

- Works Fast up to 30,000 Hz.
- Extremely quiet
- Handpiece stays cool during entire procedure
- Low volume of water is needed, making it easier to see
- Linear movement of tip – No hammer effect
- Includes three tips

Replacement Tips

- DTP00784 – Universal Tip #39
- DTP00785 – Flat Tip #32
- DTP00791 – Interproximal #37
- DTP00786 – Special / Perio Tip (optional) #33

DTP00765 – Pressurized water bottle for use with all Scalers

See page 20 for tip diagram for the
AMDENT and MAXISONIC SCALERS

TWO YEAR WARRANTY

Made in Sweden

Fine 25k Fine 30k Wide 25k



DTP10880 DTP10881 DTP10879

SWIFT "PIEZO ELECTRIC" ULTRASONIC SCALER



DTP00769

- Works fast up to 30,000 Hz.
- Low volume of water needed, only (20ml. per minute)
- Handpiece stays cool during entire procedure
- Extremely quiet
- Controls on front panel make it very easy to use
- Linear movement of tip - No hammer effect
- Comes with 3 tips (Universal - Flat - Interproximal)
- Safe to use below the gumline

- DTP00784 – Universal Tip #39
- DTP00785 – Flat Tip #32
- DTP00791 – Interproximal #37
- DTP00786 – Special / Perio Tip (optional) #33

TWO YEAR WARRANTY
Made in USA



DTP00765
Swift w/ 500 ml
Reservoir Btl.

MAXISONIC II PIEZO SCALER POLISHER



DTP00789

TWO YEAR WARRANTY
Made in Sweden

This unit allows you to scale and polish and section the teeth. It comes with two tips for scaling (one universal and one flat). It also comes with disposable prophy angles and prophy paste.

Replacement Tips

- DTP00784 – Universal Tip #39
- DTP00785 – Flat Tip #32
- DTP00786 – Special / Perio Tip (optional) #33
- DTP00791 – Interproximal (optional) #37



MAXI-SONIC III WITH RESERVOIR
DTP00766



MAXI-SONIC JET COMBO PIEZO/PROPHY JET
DTP00767



BIO-JATO BICARBONATE PROPHY
DTP00768

RAM 1000 MICRO MOTOR

COMPLETE PROPHY & RESTORATIVE UNIT
DTP103RM



With this unit you can section teeth, polish teeth, and perform restorative work. It has micro-torque control circuitry design that provides smooth constant power at all speeds. It has 5 sealed ball bearings (3 in front nose cone). It can be operated by hand or foot and goes from a crawl to 35,000 RPM. Also will go forward and in reverse. Comes complete with a prophy angle, prophy paste and cups, a contra angle with burs, and a diamond disc cutting wheel.

ONE YEAR WARRANTY
Made in USA

RAM MICRO-CAD RECHARGEABLE POLISHER
DTP104RM



ROTEX 782E MICRO MOTOR

DTP00780



Features:

- 30,000 RPM Slowspeed Handpiece
- Prophy Angle with Prophy Paste & Cups
- Contra Angles with Cutting Burs
- Variable Speed Foot Control

With these two Micro Motor units you will be able to Section Teeth, make preparations for Endo procedures, polish teeth and perform restorative work.

ONE YEAR WARRANTY

DENTALAIRE X-RAY PRODUCTS

BIO-RAY SDX

BioRay Intraoral Digital Radiography System

offers you unbeatable high quality resolution. Thicker sensor is more durable and of greater quality. 1.5mm fiber optic plate drastically improves inconsistencies and bisecting angle technique. **The software for the BioRay is made by Apteryx, a leader in software for imaging systems with guaranteed quality.** This software is very easy to use and easy to install. Installation is done by remote control of your computer (via an internet connection), which allows Sigma Biomedics to quickly set up your initial installation and troubleshoot problems if necessary. **Free updates** to the software for as long as you have the system and it bridges with virtually all practice management software. This system is upgradeable for multi-users.

Three year warranty.

Exceptional customer support and reliability at an affordable price.

Bio-Ray Intraoral Digital Radiography System Includes:

- Bio-Ray Sensor
- Bio-Ray USB Interface
- Xvlite Single User Veterinary Software by Apteryx
- Sensor Sheaths

Intraoral Sensor Specifications:

Sensor Size	Size 2
Active Pixels in Sensor	2.04M
Pixel Size (µm)	20x20
Active Image Area (mm)	34x24
External Dimensions (mm)	42.5x30.5



3 YEAR
WARRANTY

Bio-Ray System Advantages are:

- Superior software by Apteryx, easy to use, very stable
- Easy installation of software, set up "done online"
- Free updates to the software
- Built in MS Word templates for impressive client reports
- Advanced technology provides more consistency in taking x-rays.

True DICOM Compliant System Now Available

Call Dentalaire Sales for Complete details.

System Requirements

While Sigma Digital X-Ray computer requirements are actually rather modest, Sigma Digital X-Ray suggests getting the latest and best equipment you can afford for optimum performance. A higher quality screen/monitor can improve image quality substantially.

Ideal System Recommendations

Sigma Digital X-Ray suggests the following for optimum performance:

- 1.0 GHz or higher processor (faster is always better, especially for servers)
- 256 MB of RAM (again, more is always better)
- 40GB Hard Drive
- CD-R drive
- 100BT or higher Local Area Network (TCP/IP preferred)
- USB 2.0

Minimum System Requirements

Although we don't suggest these specifications for today's digital dental office, Sigma Digital X-Ray has set the following as minimums to operate our products.

- 400Mhz or higher processor
- 128MB of RAM
- CD-ROM drive
- Local Area Network for networked installations (TCP/IP preferred; A 100-BASE T network is optimum)
- USB 2.0

Sigma Digital X-Ray products are designed for Windows 9x, ME, NT4+, XP and VISTA for both servers and workstations.

For Laptop users, a power card with 3 prong grounding structure is recommended or additional hardware may be required. Call Dentalaire Sales for details.

BIO - RAY



BIORAY SYSTEM WITH #2D SENSOR
ITEM # DTPBIORAY

BIO-RAY RADIOLOGY PRODUCTS



BIO-RAY Cost/Benefit Analysis

Why pay over \$8,500.00 for a digital dental radiology system?

Some of the reasons are detailed below. The benefit analysis is based on a single doctor practice seeing five dental patients a week, and taking 20 dental films per week. This is a very conservative estimate, which results in 960 dental films per year. Most veterinarians implementing improved dental care into their practices easily average double that number of films. Multi-doctor practices enjoy multiple savings.

Item	Savings
1. Cost of Film	\$960/yr
2. Cost of Chemistry	\$336/yr
3. Employee time to develop, rinse, handle, and catalog the films. This is estimated in human dentistry to be \$2-5 per film. If this seems high, consider the total costs for your employees, which is much higher than their hourly wage. We use the lowest figure.	\$1,920/yr
4. Filing system for films. Simple envelopes or dental film mounts.	\$40-300/yr
5. Fewer re-takes. Conservatively, you will need to re-take your films around 10% of the time for improper technique. These re-takes use film, chemistry, and more employee time.	\$330/yr
Total hard cost savings per doctor per year	\$3,646/yr

There are also some "soft cost savings" that are more difficult to quantitate, yet result in undeniable positive financial effects.

1. Chemical tracking and disposal costs are a normal part of human dentistry. This will inevitably become a factor for veterinarians in the future. Registration fees, compliance paperwork and hazardous chemical disposal fees are all likely, and are not required for digital systems.
2. The client report templates built into the XV-lite software provide an easy way to present the client with pictures and radiographs of their pet. The impact of these reports on client's perception of value is substantial.
3. The software can easily allow you to display all of your patient images to clients on a separate screen in the exam room. The impact for treatment acceptance is amazing. Owners will finally see and understand the pathology in their pet.
4. No clinical downtime while you wait for a film to be developed. This can save many hours per year.
5. XV-lite is installed for you, via the Internet. You save time and hassles.
6. XV-lite can be integrated with virtually all practice management programs.
7. XV lite is more than a dental radiography program. It serves as an integrated image management suite, which allows you to gather all digital information into one place, where it can be easily accessed. Digital radiographs, digital photos, MRI's, CT's, scanned radiographs, ultrasound images and clips, endoscopic pictures and clips, etc. can all be gathered in one easy-to-access area. From there, they can be viewed, printed, and shown to clients in the exam room.

DENTALAIRE X-RAY PRODUCTS

DAZZLE YOUR CLIENTS WITH THESE IMPRESSIVE REPORTS

(Utilizing the Apteryx software in the BioRay IntraOral Radiography System)

January 26, 2006

Patient Report

Moë Wilkins

Referring doctor: Dr. Jodi Kading
B-2141



ANIMAL DENTAL CARE
Troy M. Woodard DVM, AVDC
2704 N. Academy Blvd.
Colorado Springs, CO 80918
719-535-8943



Left lower side, showing the draining tract with the root fragment underneath.



Radiograph showing the root fragment, which was extracted.



Fractured tip of the right upper canine tooth.



Left upper site prior to cleaning.

Moë presented for evaluation of dental problems. He was known to have a fractured tip of his right upper canine tooth. Us for anesthesia, several radiographs were obtained of his teeth. His two upper molars were extracted due to severe bone loss from infection. The fractured right upper canine does not appear to be infected at this time, and the fractured area was smoothed and sealed with bonded resin. He was seen to have a draining tract over the area where his left lower first premolar used to be. The drainage was associated with a retained root fragment, which was extracted. The surgical site was closed with absorbable suture material. As a courtesy, his teeth were also cleaned today and several areas of root exposure were treated with closed root planing. Please care for him as follows:

1. Give one ml of the Aesirbe drops AM and PM, for 10 days, for infection.
2. Give 1 small syringe (0.1 ml) of the Suprenex solution, AM and PM, as needed for pain. This may make him drowsy, and the dose may be halved or doubled as needed.
3. The salines will dissolve over the next few weeks.
4. We would like to re-check his mouth in 30 days to assess healing.
5. Moë would benefit from any home care you are able to do.
6. Please call with any questions or concerns.



SIGMA BIOMEDICS XV-LITE SOFTWARE

DTP61400

MINI X-RAY VIEWER



DTP10500 Colors: Grey & Blue

DTP10501 Replacement Bulb

X-RAY ITEMS

- DTP10439 – KODAK X-Ray Film #0, box of 100 (DF-54)
- DTP10440 – WOLF X-Ray Film #2, box of 100 (D-Speed)
- DTP10441 – KODAK X-Ray Film #2, box of 150 (DF-58)
- DTP10455 – KODAK X-Ray Film #4, box of 25 (DF-50)
- DTP10431 – Film Carriers (for film #2) Pk. of 25
- DTP10432 – Film Carriers (for film #4) Pk. of 25
- DTP10460 – X-Ray Film Clips, Pk. of 3 (Useful in film developing)
- DTP10430 – X-Ray Film Holder (autoclavable)
- DTP10433 – X-Ray Film Mounts #2 (100) (For film storage)
- DTP10434 – X-Ray Film Mounts #4 (25) (For film storage)
- DTP10484 – X-Ray Developer (quart)
- DTP10486 – X-Ray Fixer (quart)
- DTP10485 – X-Ray Fixer/Developer (set of 2)
- DTP10442 – Size 1 X-Ray Film (100) KODAK
- DTP10443 – Size 3 X-Ray Film (100) KODAK

CHAIRSIDE DARKROOM



Complete with:

- X-Ray Fixer
- Developer
- Film clips (Set of 3)
- Convenient
- Easy to use
- Cost effective

DTP10115

DTP10114 –
Replacement Jars
(4) for Chairside
Darkroom.

DTP10113 –
Replacement Red
Cover for Chairside
Darkroom.

DENTALAIRE X-RAY PRODUCTS

Corix® Pro 70
Vet Dental X-Ray

With the most powerful and user friendly hand held controller in the market that fits into the palm of your hand.



Corix Pro 70 Wall Mount, DTP60800

• All is needed is to select the tooth and patient size. A powerful microprocessor will automatically choose the exposure time, according to the type and position of the tooth, patient size, film or digital sensor, with automatic compensation for the main voltage variation, or manual override if required. Efficiently obtain the optimized technical factors that will provide constant radiographic density, with less radiation exposure. Full protection against x-ray tube overheating: the microprocessor will set the minimum cooling period between exposures, thus extending the tubehead life.

• **Standard, Short or Long Extension Arm,** to fit all your requirements.

• **Long reach Double Positioning Arm**
Moves effortlessly and folds compactly when not in use. Ease of positioning and drift free stability.

• **70 KVp Tubehead**
Assures high quality, sharp, well defined images in intra-oral radiography.



• **Corix® Pro 70**

Wall mount, also for side-by-side operatory stations.

Technology and reliability in X-Ray equipments, since 1974

• **Corix® Pro 70**

Mobile Stand, on a sturdy and stable rolling stand.

Corix Pro 70 Mobile Stand, DTP61300



DENTALAIRE X-RAY PRODUCTS


70 KVp Tubehead

The 70KVp-8mA Tubehead with a small focal spot, low energy filter and beam limiting device assures outstanding sharpness and maximum details in intra-oral radiography. Positioning is made easier with an angle indicator on both sides of the unit. Lead shields and critical components are oversized to provide more years of safe and trouble free operation.



Hand Held Digital Exposure Timer

A new generation of small size, 3 digit remote controller with a 1/100 of a second resolution, fully compatible with the fastest film or a filmless digital imaging system. 9 factory pre-set memories, patient size, film or digital sensor, or manual override if required. Automatic compensation for the main voltage variations. All of above and more in a controller that fits into the palm of your hand. Optional Remote Timer Configuration, should the controller be located outside the operatory.

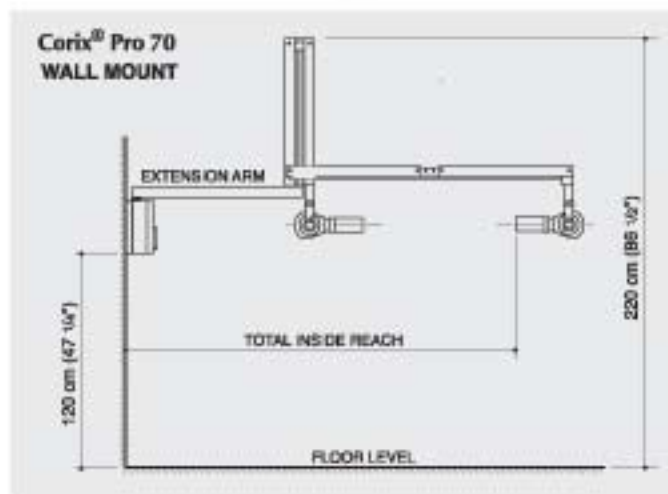

Long Reach Double Positioning Arm

The long reach scissors arm moves effortlessly through a wide area, yet folds compactly when not in use. Ease of positioning and drift free stability.

Standard, Short or Long Extension Arm (Wall Mount version)

Designed to fit your requirements:

Standard Extension Arm provides a total inside reach of 180.6 cm. (71 1/8")

Short Extension Arm provides a total inside reach of 135.6 cm. (53 3/8")

Long Extension Arm provides a total inside reach of 190.5 cm. (75")

Corix[®] Pro 70 is a wall mounted or mobile system

Conditions in your practice will dictate where to perform x-ray examinations. With the Corix Pro 70 you can choose the option that best fits your needs:

Corix[®] Pro 70 Wall Mount. DTP60800

Wall mounted in the x-ray room or directly in the treatment room, even for side-by-side operatory stations.

Corix[®] Pro 70 Mobile Stand. DTP61300

Mobile stand, to reach separate operatory stations on a sturdy and stable rolling stand.

1 Year Warranty. Price does not include installation.

Technical Features

RATED LINE VOLTAGE	120V +/- 10% (108-132 V a.c.), 50/60 Hz.
POWER CONSUMPTION	1.050 KW max @ 132 V a.c.
RATED TUBE POTENTIAL	70 Kv +/- 7% (Single phase, self rectifying)
RATED TUBE CURRENT	8 mA +/- 15% @ 120 V a.c.
EXPOSURE TIME RANGE	0.03 to 3.00 s., in steps of 0.01 s. Plus 9 pre-set exposure times, with automatic main voltage compensation. A pre-heating time of 0.15 s. for loading the x-ray tube is not shown on the display.
MINIMUM FOCAL SPOT TO SKIN DISTANCE	200 mm. (8")
RADIATION LEAKAGE @ 1 m.	Less than 50 mR/h (Technical Factors: 70KVp, 8 mA, 1 s.)
FOCAL SPOT	0.6 x 0.6 mm.
TOTAL FILTRATION	2 mm. Al eq.
NET WEIGHT	Mod. WM (Wall Mount): 28 Kg. Mod. MM (Mobile Stand): 40 Kg.
SHIPPING WEIGHT	Mod. WM (Wall Mount): 33 Kg. Mod. MM (Mobile Stand): 46 Kg.
SHIPPING DIMENSIONS	Mod. WM (Wall Mount): 1 Box, 131 x 52 x 26 cm. Mod. MM (Mobile Stand): 2 Boxes, 111 x 34 x 19 cm., and: 86 x 34 x 33 cm.
<i>Corix Medical Systems[®] reserves the right to modify the design and specifications without prior notice.</i>	

[®]Corix, and: Corix Medical Systems, are Registered Trademarks.

Each CORIX comes with a book on veterinary dental radiography, as well as technique chart

DENTALAIRE X-RAY PRODUCTS

LEAD APRON



.50mm (Middleweight),
24" x 36", Many colors available.
Most popular colors: Burgundy,
Teal, Royal Blue
DTP63001

.50mm (Middleweight),
24" x 42", Same colors available
as above apron.
DTP63011

THYROID COLLAR



.50mm provides protection to the neck
area. The collar is in place by a Velcro
closure. Most popular colors: Burgundy,
Teal, Royal Blue.
DTP75060

COMFORT FLEX X-RAY GLOVES

.50mm, a five finger molded glove which offers even more
flexibility than normal gloves. Excellent for any application
requiring dexterity for small detailed work. The exterior finish
can be cleaned with soap and water. Sizes available: Small,
Medium, Large.
DTP12406 (pair)



- More flexibility than normal gloves
- Excellent for any application requiring small detailed work
- Exterior finish is easily cleaned with soap and water



HANDPIECE PARTS & ACCESSORIES

PUSH BUTTON CONVERSION CARTRIDGE

Changing burs is as easy as pushing a button!



DA400RC – Standard
Replacement Cartridge
DA400CP – Cap Only



DA400RCA – Push button
replacement turbine with cap
DANSK-A485 – Turbine Only



DA310BT – Handpiece
Bur Changer



DA400FOB – Fiberoptic Bulb



DA400CH – High Speed
Replacement Chuck for
standard turbine



DA2500SNC - E-Type
Straight-nose cone
20,000 RPM



DA2500ECA - E-Type
Contra-Angle, Latch type
20,000 RPM



DA2500EM - E-Type Motor
20,000 RPM

HIGH SPEED HANDPIECES



DA400A – Standard High Speed Handpiece



DA400APB – High Speed Handpiece with push button



DA400AFOPB – Push Button High Speed Fiber Optic Swivel Handpiece



DTP400AFO – Fiberoptic High Speed Handpiece

LOW SPEED HANDPIECES



DA600HA – Low Speed Handpiece Hygienist (5000 RPM) with push/pull chuck system



DA600A Light Weight - Low Speed Handpiece Standard (5000 RPM) with opening & closing chuck system and reverse



DA557FG – Contra Angle Attachment



DA553UL – Contra Angle Attachment



DA700A – Low Speed, Light Weight, Lube Free Handpiece

SCALER HANDPIECE



SONIC AIR SCALER HANDPIECE
DA6000AS – Sonic Air Scaler with 3 Tips



DA100PL – PEN-LUBE Handpiece lubricant with fine needle tip applicator 1/4 oz.

DA400GAS – Gasket for 4-hole handpieces

DA500GAS – Gasket for fiberoptic handpiece

DA400FOB – Fiber Optic Bulb Replacement



DA100SL – ONCE-A-DAY Handpiece Lubricant 8 oz. can

VETERINARY DENTAL MODELS





For Dental Education and Demonstration

Transparent and firm smooth plastic models with removable and replaceable teeth. Graphically shows important anatomical relationships.

DTP101VMD Canine DTP101VMDCO Canine w/Gingiva
 DTP102VMC Feline DTP102VMCCO Feline w/Gingiva
 DTP103VMDC Combo Canine and Feline

DENTAL CHART STICK-ON LABELS

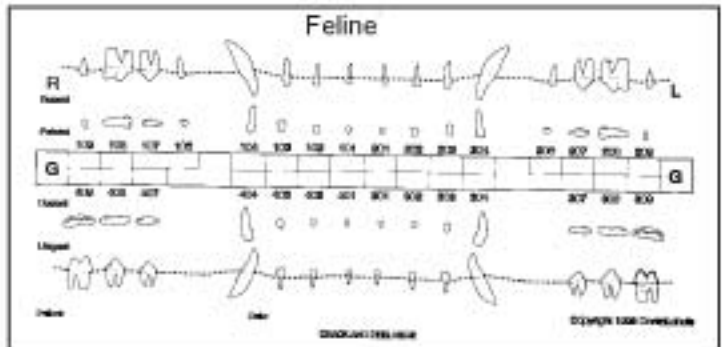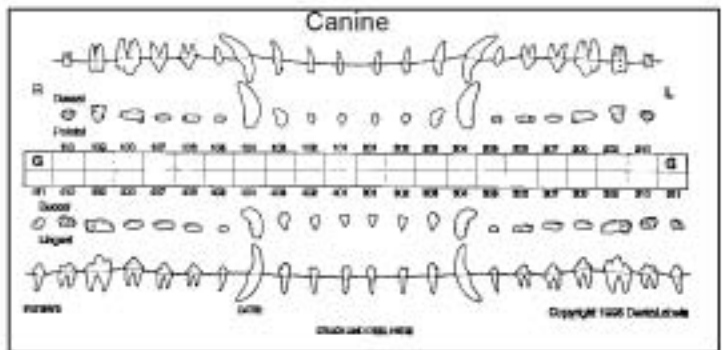

DTP20340 – Starter Kit - 25 Canine / 25 Feline
 DTP20341 – Canine Refill - 50
 DTP20342 – Feline Refill - 50

SHARPENING SYSTEM

The only sharpening system designed to maintain and restore instruments to a like new condition. Save time and \$\$\$ on your surgeries. Sharp instruments make your job easier. This machine will properly sharpen: Elevators, Curettes, Scissors, Olsen Hegar Needle Holders, Orthopedic Instruments and much more. Most instruments can be sharpened in 2 minutes or less. Used in dental practices, hospitals, and in veterinary clinics for over 30 years. You can't afford not to have it. Give it a try! You'll be glad you did.



DTP90602–
 RX Sharpening
 System II

3 YEAR
 WARRANTY

Made in USA



RX HONING WORK STATION

Includes Lamp,
 Magnification and Base

DTP90601

OptiVISOR



DTP90608 – DA-2
 (1-1/2 x at 20")

DTP90609 – DA-3
 (1-3/4 x at 14")

DTP90610 – DA-4
 (2 x at 10")

DTP90612 – Light
 for visor

OPTICAL GLASS BINOCULAR MAGNIFIER

- Leaves both hands free
- Optical polished (not plastic) glass lens
- Tilts up when not in use
- Adjusts to all head sizes • Leather comfort band
- Can be worn over prescription or safety eyeglasses

RX HONING FLEX SHAFT

Heavy duty shaft
 provides pencil
 control in the honing
 of intricately shaped
 cutting edges.



VETERINARY DENTAL WET LABS

Even the best equipment requires some knowledge to use effectively. Dentalaire is pleased to recommend the dental training classes presented by Dr. Tony M. Woodward DVM, AVDC. Dr. Woodward is board certified in veterinary dentistry, and has trained hundreds of doctors and technicians in practical veterinary dental techniques. His experience in starting multiple veterinary practices makes him uniquely qualified to teach those techniques most useful for small animal practitioners. His facility is located one hour south of Denver, in beautiful Colorado Springs, Colorado. High instructor to

student ratios and excellent equipment ensure a quality training experience for all participants. In just two days, you can learn all of the techniques needed to deliver 90% of all the dental problems encountered in small animal practice.

Eight CE credits are available for each day. Classes are a combination of lecture and wet lab, and include handouts, breakfast, lunch, lab materials, individual specimens, and complimentary instrument sets. These two classes cover the skills required to deliver 90% of the dental care required in small animal practice.



Course Content:

Level I is for both technicians and doctors, is limited to 14 participants, and covers dental cleaning, oral pathology, instrumentation, equipment, dental radiology, and digital dental radiography systems. This class will give you the skills required to provide a good cleaning procedure and identify oral pathology. The wet lab portion focuses on obtaining diagnostic dental radiographs.

Level II is for doctors only, is limited to seven participants, and covers surgical extraction technique, introduction to periodontal surgery, oronasal fistula repair, and bonded sealants for fractured teeth. Three instructors for seven participants ensure a high level of individual attention. Attend both classes and save \$25.00



Classes are held on Friday and Saturday, with a level one class followed by a level two class the next day. This allows interested individuals to attend both classes in one weekend. Classes are scheduled throughout the year.

For more information and class availability, contact Animal Dental Care at:

(719) 536-9949

Experience Excellence

The Ultimate Veterinary Dental Station

The only unit you will need for all of your dental procedures... With the Dentalaire "Ultimate" Dental Station there are no "options" needed. The unit is totally complete with everything you need to perform all of your dental procedures.


Features:

- Virtually soundless compressor
- Two handpiece control system
- Tank powder coated inside to prevent corrosion
- Dual filtration system, which gives you cleaner air
- Three way (air, water and spray) syringe
- Built in dual water system with (2) 2 liter water bottles and inline water filter
- Built-in ultrasonic scaler with autoclave tips
- High speed fiberoptic handpiece with push button turbine
- Low speed polishing handpiece with prophylaxis angle prophylaxis paste and cups
- Telescopic height adjustment from 38" down to 26" for convenient storage
- Work area 22" x 18" with stainless steel removable tray
- Larger 9 outlet multiple strip for accessories



High Speed Fiberoptic Handpiece

with oral illumination system. Push button control for ease changing.

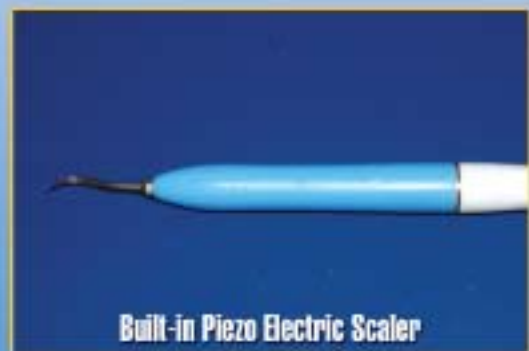

Built-in Piezo Electric Scaler

with Universal and Flat Scaling tips. Runs cool and quiet at 30,000 cycles per second.



Light Weight, Low Speed, Lube Free Handpiece

with "reciprocating" prophylaxis angles. No more heat - no more splatter and no more tangles.



800.866.6881 • www.dentalaireproducts.com

Meeting and Exceeding the Needs of
Veterinary Dentistry Since 1988