


# DENTALAIRE X-RAY PRODUCTS

## Corix® Pro 70 Vet Dental X-Ray

*With the most powerful and user friendly hand held controller in the market that fits into the palm of your hand.*

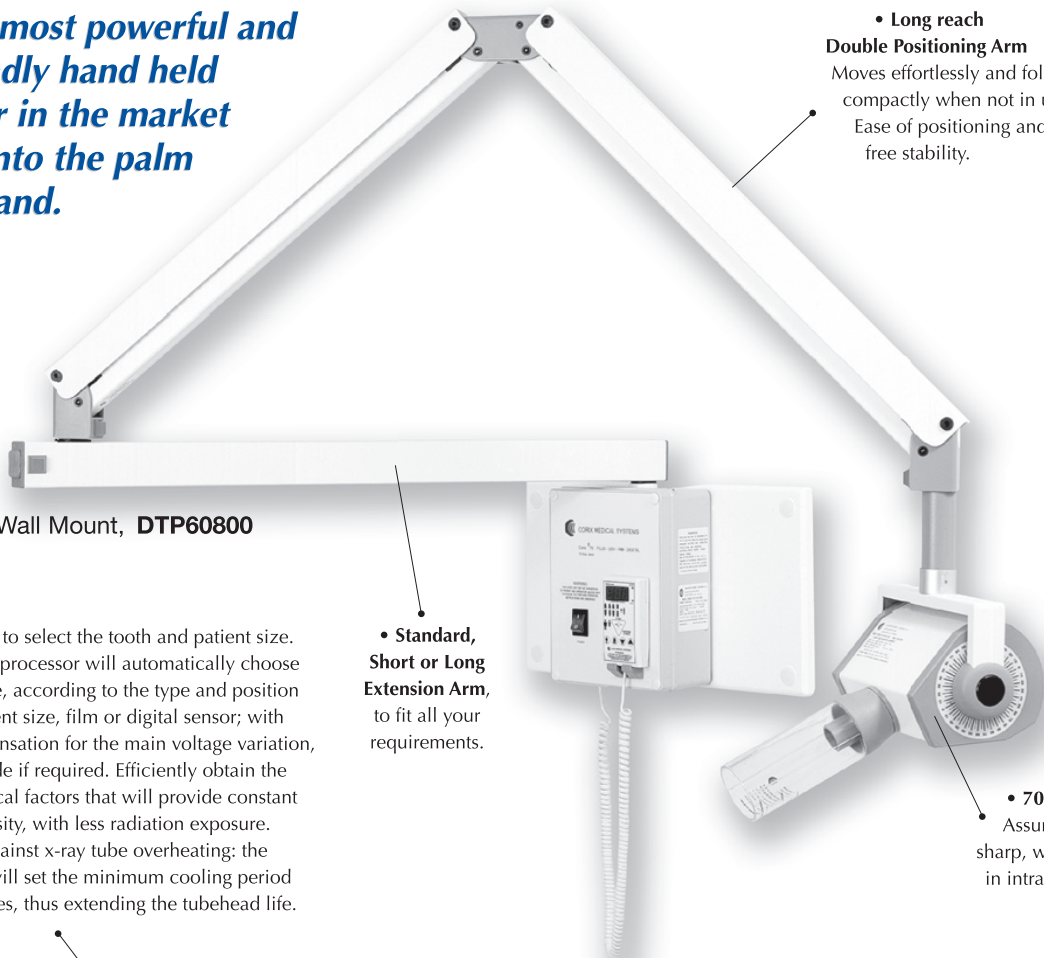

Corix Pro 70 Wall Mount, **DTP60800**

• All is needed is to select the tooth and patient size. A powerful microprocessor will automatically choose the exposure time, according to the type and position of the tooth, patient size, film or digital sensor; with automatic compensation for the main voltage variation, or manual override if required. Efficiently obtain the optimized technical factors that will provide constant radiographic density, with less radiation exposure. Full protection against x-ray tube overheating: the microprocessor will set the minimum cooling period between exposures, thus extending the tubehead life.

• **Standard, Short or Long Extension Arm**, to fit all your requirements.

• **Long reach Double Positioning Arm**  
Moves effortlessly and folds compactly when not in use.  
Ease of positioning and drift free stability.

• **70 KVp Tubehead**  
Assures high quality, sharp, well defined images in intra-oral radiography.

### • **Corix® Pro 70**

*Wall mount, also for side-by-side operatory stations.*

**Technology and reliability in X-Ray equipments, since 1974**

### • **Corix® Pro 70**

*Mobile Stand, on a sturdy and stable rolling stand.*

Corix Pro 70 Mobile Stand, **DTP61300**

