

Meeting and Exceeding the Needs of Veterinary Dentistry since 1988 www.DentalaireProducts.com 800.866.6881

The Ultimate Dental Station

The only unit you will need for all of your dental procedures!

With the Dentalaire "Ultimate" Dental Station there are no "options" needed. The unit is totally complete with everything you need to perform all of your dental procedures.

Features:

Virtually soundless compressor

Oiless SilentAire compressor available

Two handpiece control system

Tank powder coated inside to prevent corrosion

Dual filtration system, which gives you cleaner air

Three way (air, water and spray) syringe

Built in dual water system with (2) 2 liter water bottles and inline water filter

Built-in ultrasonic scaler with autoclave tips


High speed fiberoptic handpiece with push button turbine

Low speed polishing handpiece with prophy angle prophy paste and cups

Telescopic height adjustment from 38" down to 26" for convenient storage

Work area 22" x 18" with stainless steel removable tray

Larger 9 outlet multiple strip for accessories

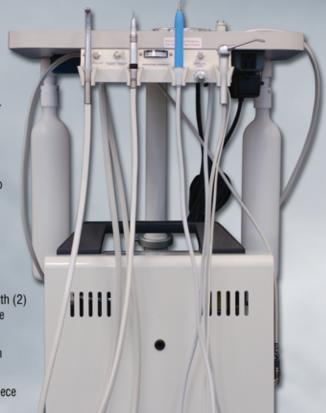

DENTAL MIRES



High Speed Fiberoptic Handpiece

with oral illumination system. Push button control for ease changing.



Built-in Piezo Electric Scaler

with Universal and Flat Scaling tips.
Runs cool and quiet at
30, 000 cycles per second.



Light Weight, Low Speed, Luie Free Handpiece

with "reciprocating" prophy angles. No more heat - no more splatter and no more tangles.



Ultimate Dental Station.