

Dentalaire Prestige Dental Station

Mobile Dental Station



Prestige Dental Station

The **Prestige Dental Station** is our most popular mobile dental station. It is known for its reliability, internal storage, and pneumatic height adjustment.

The included handpieces offer versatility and efficiency. This dental station now includes the **NEW** Swivel LED High-speed Handpiece that provides a brighter light source, compared to that of fiber-optic tubing.

- Quick disconnect pressurized water connection
- NEW High-speed LED swivel handpiece
- Low-speed handpiece
- LED Piezo scaler
- Telescopic pneumatic height adjustment from 38 1/4" - 30"
- 2 - year warranty excluding handpieces and accessories
- All other warranties are described in the individual product description

DTP00511 Prestige Dental Station with Oil Compressor
DTP00511-OL Prestige Dental Station with Oil-less Compressor



Internal Storage and Parts

Opening the lid, will give access to the internal storage and parts.

- Easy access to replaceable parts
- Storage tray
- Improved serviceability



High-speed LED Swivel Handpiece

The High-speed Swivel LED Handpiece provides direct light on the face of the handpiece.

- Eliminates fiber optic tubing
- Smooth and quiet operation



Adjustable Height

With a pneumatic lift, the height for this dental station can be easily adjusted.

- Dimensions: 38 1/4" - 30" H x 22 1/2" W x 24 1/4" - 22" D
- Easy to store when not in use


Innovative design.

Included with the Prestige Dental Stations is the Dentalaire Vet Pack, which is a great starter kit for using your dental station. This kit includes:

- Mint prophy paste
- Snap-type prophy cups
- Snap-on prophy angle eco.
- Bur pack
- O-ring for tank drain
- All-In-One handpiece cleaner and lubricant



Looking to customize your dental station?
Contact us to learn more about your options.

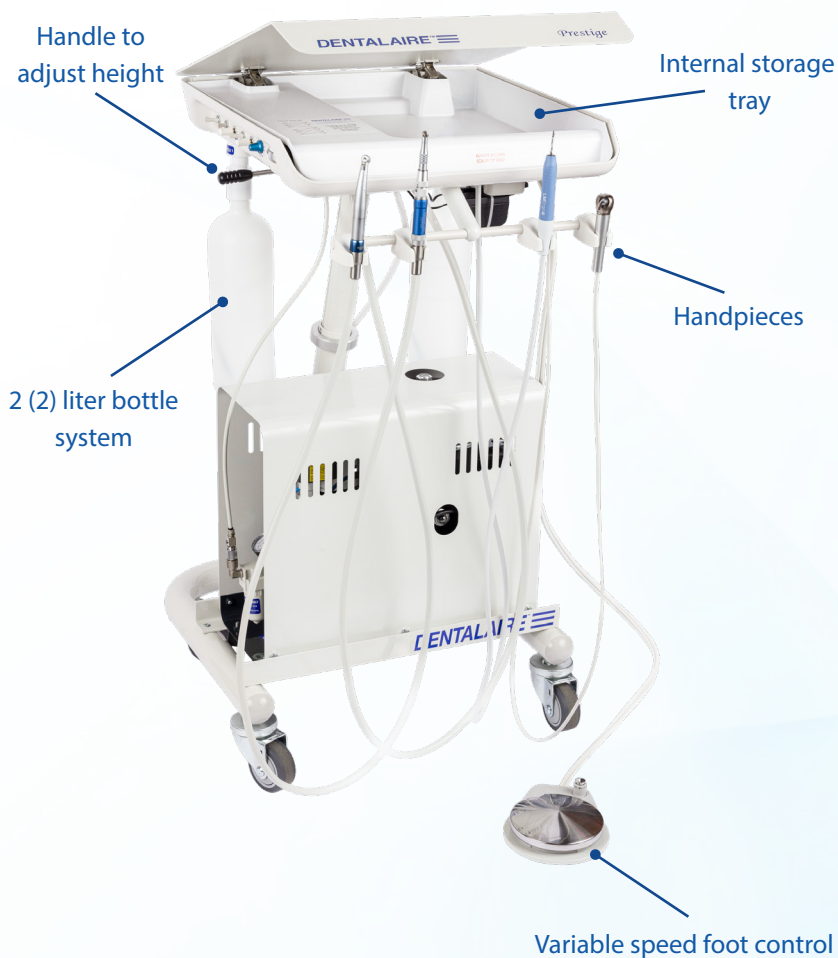

Oil Compressor

- Horse Power: 1/2
- Noise Level: 40 decibels
- Tank Size: 2.25 gallons



High-quality Handpieces

- High-speed LED Swivel
- Low-speed
- LED Ultrasonic Scaler



Side Panel Controls

- Asepsis design
- User-friendly controls
- Smooth handpiece tubing